





Poli VAC (Polish Device)

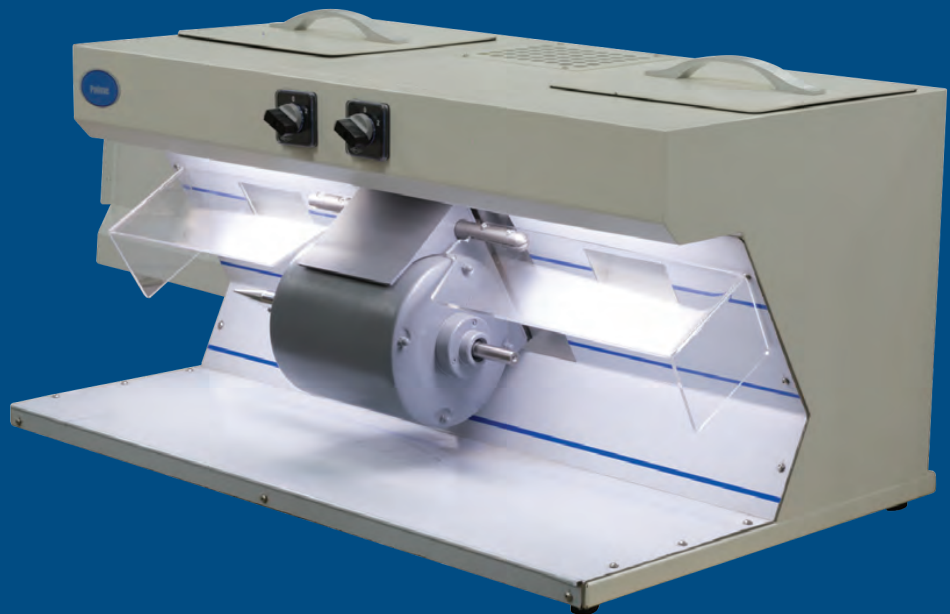

QUALITY & ADVANCED TECHNOLOGY
PROVIDE FULL FACILITIES
AS YOU DESIRED

Poli VAC (Polish Device)

Polish Device with high strength for polishing & processing the surface, along with strong suction results in a high quality disbursement, along with maximal hygiene. So the user will reach to maximum efficiency of this product

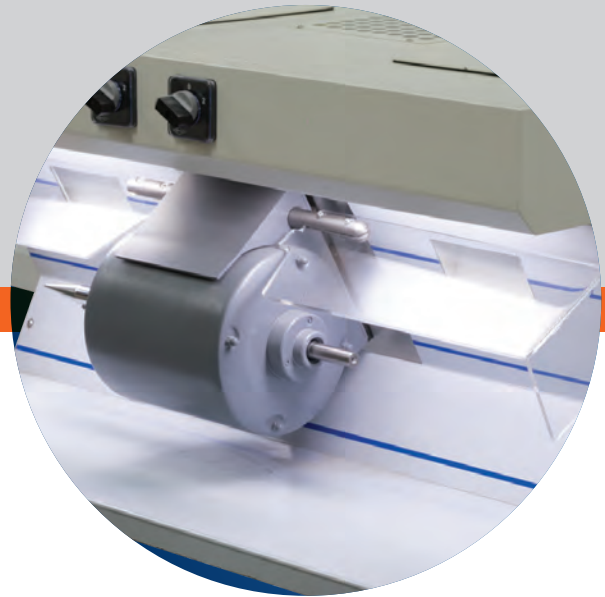

Features

- Usable for all removable prostheses: acrylic and chrome-cobalt dentures
- Final processing ability with different rotation
- Two installment axes for processing and polishing brushes
- Automatic suction system in order to prevent suspended dust
- Fluorescent light system for better visibility
- Suction bag for easy evacuation
- Plastic reservoirs for easy cleaning
- Mechanical system for extracting brush from shaft for ease of performance & maintenance

Technical Data

- Dimensions: 51 cm (width)×81 cm (depth)×45 cm (height)
- Net weight: 55 Kg
- Maximum power input: 500 W
- Operating voltage: 210-240 V-50/60 Hz
- Motor Angular velocity: 400-1200 RPM

Address: Humboldt Street 12 - D40723
Hilden - Germany
Tel: + 49 (0) 21033373500
Fax: + 49 (0) 21033373550
Email: info@polarisgmbh.de
www.polarisgmbh.de

24
Months Warranty

10
Years Supporting of
Spare Parts